

CTD- χ pod Processing Guide

Andy Pickering

September 16, 2017

Contents

1	Introduction	3
2	Software	3
3	Data Organization	3
3.1	Raw Data	3
3.2	Processed Data	4
4	CTD data	4
4.1	24Hz CTD Data	4
4.2	1-m Binned CTD Data	4
5	Brief Outline of Major Processing Steps	5
6	Processing steps	6
7	χpod Processing Parameters	7
8	Output File Formats	8
8.1	/cast	8
8.2	/cal	8
8.3	/avg	9
9	Notes	10
10	Post-processing Analysis	10
11	Sensitivity Analysis	10
12	Standard Plots made by processing	10

13 Example Figures	12
14 m-files	21

1 Introduction

This is a guide to processing data from CTD- χ Pods. χ Pods are self-contained instruments developed by the OSU Ocean Mixing Group to measure turbulence. They have a FP07 thermistor that measures temperature and temperature gradient, sampled at 100Hz.

2 Software

- The processing is done with Matlab codes, which are maintained in github repositories:
- <https://github.com/OceanMixingGroup/mixingsoftware>.
- The CTD- χ Pod specific codes are in `/mixingsoftware/CTD_Chipod`.
- Codes for processing raw CTD data (.hex format) are contained in `/mixingsoftware/ctd_processing`.
- Code for processing specific cruises is kept in the 'cruises' repository:
https://github.com/OceanMixingGroup/cruises/tree/master/ctd_chipod

See 'Processing Steps' section for more details. Before processing, you should clone the 'mixingsoftware' and 'cruises' repositories and make sure you have the latest versions of both. The processing is designed so that code, figures and notes are kept in the 'cruises' repository. The data are too large to store on github, and must be stored on a server or local machine. A script in the processing specifies this location.

3 Data Organization

Raw data should be organized into folders as listed below. It is important to use this standard structure so that the processing functions can be applied to each cruise. The function `make_ctd_chipod_folders.m` will make the empty folders when you start processing a project. `Project` is the name of the cruise ('IO8' for example).

3.1 Raw Data

- `[Project]/Data/raw/Chipod/SNxxxx` : Raw χ pod data.
- `[Project]/Data/raw/CTD/` : Raw CTD data (.hex,.XMLCON,.hdr).

Processed data is organized as follows:

3.2 Processed Data

- [Project]/Data/proc/Chipod/ Main folder for processed χ pod data
- Data for each χ pod is stored in /[whSN]/[pathstr] where **pathstr** is constructed based on Params, for example:
`zsm20m_fmax7Hz_respcorr0_fc_99hz_gamma20`
This allows us to run the processing for different Params and easily compare the results.
- [Project]//Data/proc/CTD Processed CTD data. (/24Hz,/binned).

4 CTD data

CTD data needs to be obtained and processed before running the χ pod processing. Processing requires:

1. 24Hz data: used to determine the time offset between the CTD data and χ pod data, by aligning (highpassed) dp/dt from the CTD with w estimated by integrating vertical chipod accelerations. 24Hz temperature is also used to calibrate the chipod temperature voltages and convert to physical units.
2. 1-m binned data: used to calibrate chi-pod T, and calculate dT/dz and N^2 .

I prefer to use codes in /mixingsoftware/ctd_processing/ to process standard Seabird hex files into standard mat files to use in chi-pod processing. See `Process_CTD_hex_Template.m` for an example. The data can be processed using other methods, but must contain the specific fields required by the processing routines as listed below:

4.1 24Hz CTD Data

Required fields for 24Hz CTD data are:

- dnum : datenum
- t : temperature [C]
- p : pressure [db]

4.2 1-m Binned CTD Data

Required fields for binned CTD data

- dnum : datenum

- t : temperature [C]
- p : pressure [db]
- s : salinity [psu]
- lat : Latitude [$^{\circ}$]
- lon : Longitude [$^{\circ}$]

5 Brief Outline of Major Processing Steps

1. Process raw CTD casts if necessary
2. Find chipod data for each CTD cast
3. Align chipod data with CTD data
4. Calibrate chipod data
5. Apply chipod calc in 1 sec windows

More details are given in the next section.

6 Processing steps

1. Clone the ‘mixingsoftware’ and ‘cruises’ repositories and make sure you have the latest versions of both. If you have cloned them already, you update to the latest version with `+git pull+` from terminal, or the ‘sync’ button in Github Desktop.
2. In the `process_[project]_script`, specify the path to your ‘mixingsoftware’ folder (`mixpath`), and make sure `the_project` is correct.
3. Make/modify `Load_chipod_paths_[Project].m`. This specifies filepaths to raw data, as well as where output will be saved. You should only need to specify the base directories (assuming data is in the standard structure):
 - (a) `BaseDir` Path to directory where (raw) data is stored. Recommend making a copy of the raw data to process from, to avoid altering or deleting any of the original raw data.
 - (b) `cruise_path` Path to your local clone of the ‘cruises’ repository.

Figures, notes, summaries etc. will be saved to folders under `cruise_path`. Processed data will be saved to folders under `BaseDir`.

4. Make/modify `Chipod_Deploy_Info_[Project].m` with info for specific deployment. Much of the rest of the processing uses information from the above two files. *The info for each unit can be found in the excel log file*.
5. Plot the raw chi-pod data (`PlotChipodDataRaw_General`). This will let you see quickly if any sensors or files are obviously bad. If there are bad files you can specify these in a `Bad_File` list so they will not be loaded during processing (this is not necessary, but will reduce time loading bad files, and may prevent bad data from accidentally being used).
6. Run `MakeCasts_CTDchipod_function.m` for a few casts. Check time-offset / alignment (figure 4); if not right, probably need to modify `az_correction` or `flip_ax_az` in `Chipod_Deploy_Info`. * Note if χ pod clock is off by a large amount (for example you forgot to set clock), processing will not work because it finds the files for each cast by the timestamp in their name. You will need to apply a time offset to the raw files (make copies first) and rename them in this case.
7. Run `MakeCasts_CTDchipod` for one good cast for each SN. Check time-offset / alignment (figure 4); if not right, probably need to modify `az_correction` or `flip_ax_az` in `Chipod_Deploy_Info` script. * Note if χ pod clock is off by a large amount (for example you forgot to set clock), processing will not work because it finds the files for each cast by the timestamp in their name. You will need to apply a time offset to the raw files (make copies first) and rename them in this case.

8. Once that is figured out, you can run `MakeCasts_CTDchipod` for all casts. This script finds χ pod data for each cast, aligns the data, calibrates the data, and saves a file with data for each cast. Data are saved to
`/Data/proc/Chipod/cast` and
`/Data/proc/Chipod/cal`.
 If Matlab crashes before finishing (which seems to happen relatively often on my laptop, you can load `proc_info.mat` and look at the `last_icast` field to see where to restart `MakeCasts_CTDchipod`.
9. `Plot_TP_profiles_EachCast_CTDchipod`
10. Run `SummarizeProc`. This makes some summary tables and figures from the ‘Xproc.mat’ data saved during `MakeCasts_CTDchipod`. The two latex tables it produces can be copied and pasted into the Latex notes template.
11. `Plot_TP_profiles_EachCast_Template.m` For each cast, this plots the temperature derivative from all chi pod sensors on the same scale so you can compare them and check for bad profiles. Figures are save to `/Figures/TPprofiles/`. This is another good place to check that the up/down designations are correct; the title for the designated direction will be green. The corresponding profile should be less noisy than the other direction.
12. Run `do_chi_calc_ctd_chipod.m`. This does the chi calculation for all the cast files made in `MakeCasts_CTDchipod`. Data are saved to `/Data/proc/Chipod/avg/`
13. Modify and run `make_combined_chi_struct`. This will combine all the processed ‘avg’ files for each cast into a single structure ‘XC’.
14. Modify and run `plot_XC_summaries` to make summary plots from the combined ‘XC’ structure.
15. `Write_Latex_Notes` Makes an automated latex notes file for this cruise.

7 χ pod Processing Parameters

Variable parameters for the processing are specified in a ‘Params’ structure.

- `fmax` : Maximum frequency to integrate temperature gradient spectrum to (default 7 Hz). In practice this varies with each sensor.
- `gamma` : Assumed mixing efficiency (default 0.2).
- `resp_corr` : Option to apply frequency response correction (default 0).

- `fc` : Cutoff frequency for response correction (if applied) (default 99 for naming purposes).
- `z_smooth` : The distance [m] over which N^2 and dT/dz are smoothed (default 10m).

8 Output File Formats

8.1 /cast

‘Raw’ (ie uncalibrated) χ pod data for each cast. Structure ‘cast’ with fields:

- `datenum`
- `T1` = Temperature [V]
- `T1P` - dT/dt [V]
- `AX` - X-acceleration [V]
- `AZ` - Z-acceleration [V]
- `chi_files` - List of raw chipod files from which data was combined into this structure.
- `time_range` - Time range of CTD cast.
- `castname` - Name of CTD cast.

Ex. filename :

8.2 /cal

Calibrated data for each cast. Up and down casts are saved in separate files. Structure ‘C’ with fields:

- `datenum`
- `P` [db]
- `T1P` - dT/dt [Ks^{-1}]
- `fspd` - [ms^{-1}]
- `castidr` - up/down
- `info`
- `ctd`

Ex. filename : `rh10011_SN2013_downcast.mat`

8.3 /avg

Result of χ pod calculations. Contains data points for each 1 sec window.

- Params
- datenum
- P - [db]
- N2 - [s^{-2}]
- dTdz [Km^{-1}]
- fspd - [ms^{-1}]
- T [$^{\circ}C$]
- S [psu]
- theta
- sigma
- chi1- χ [K^2s^{-1}]
- eps1 - ϵ [Wkg^{-1}]
- KT1 - K_T [m^2s^{-1}]
- TP1var
- samplerate - [hz]
- nu
- tdif
- lat - [$^{\circ}$]
- lon - [$^{\circ}$]
- castname
- castdir
- Info

Ex. filename : avg_08402_downcast_SN1013_T1.mat

9 Notes

There is a template for standard latex notes that can be modified for each cruise/deployment: `Chipod_Notes_Template_AP`.

10 Post-processing Analysis

- Plot χ pod time-offsets. The time-offset should drift linearly with time, and may be reset during the cruise. Typically should be less than 20 sec (depends on how accurately time was set). In some cases, the time was set wrong (off by a day, wrong time zone etc.); in this case the processing script needs to be modified to correct this.
- Compare different sensors and cast directions: The quality of data obtained depends heavily on the setup of the CTD and the χ pod mounting. Typically the upward looking sensors on the upcasts are the cleanest data.
- Sensitivity to parameters: see details below.

11 Sensitivity Analysis

The sensitivity of the results to different parameters should be examined, including `fmax` and `z_smooth`.

12 Standard Plots made by processing

- `castname_Fig0_RawCTD` (Figure 1).
- `castname_Fig1_RawChipodTS` (Figure 2).
- `castname_Fig2_w_TimeOffset` (Figure 3)
- `castname_Fig3_w_TimeOffset_Zoom` (Figure 4). Shows a zoomed-in period of the time-offset. The two timeseries should match in the lower panel.
- `castname_Fig4_w_dTdtSpectraCheck` (Figure 5). This compares the analog dT/dt (T1P) to the differentiated T signal. The spectral amplitudes should agree up to a certain frequency where the digital dT/dt becomes dominated by noise. If the spectra below the noise level don't match, the time-constant used in calibration might need to be modified.
- `castname_Fig5_T_P_dTdz_fspd` (Figure 6). Summary of the aligned and calibrated data.

- `castname_Fig6_downcast_chi_T1_avgPhist` (Figure 7).
- `castnameFig7__upcast_chi_SN600_T1_avg_chi_KT_dTdz` (Figure 8). Summary of the final data ('avg' structure).

13 Example Figures

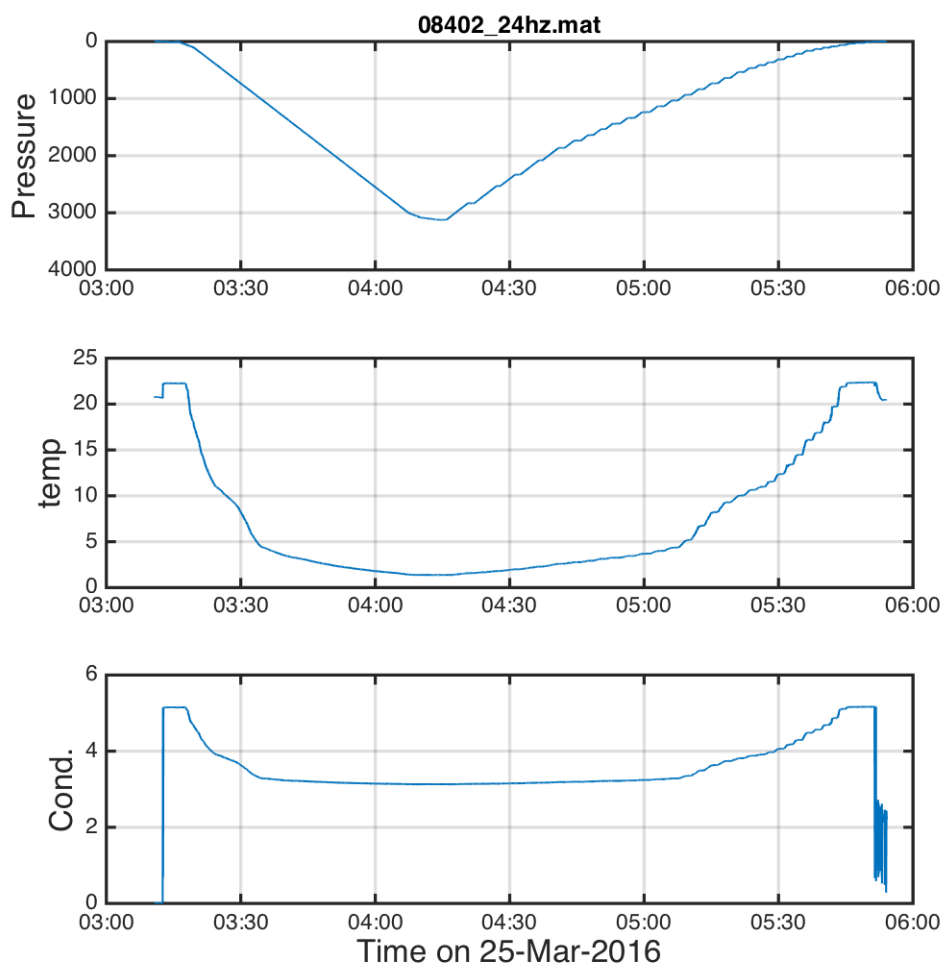


Figure 1: Time-series of raw CTD data for one cast: (a) Pressure (b) Temperature (c) Conductivity.

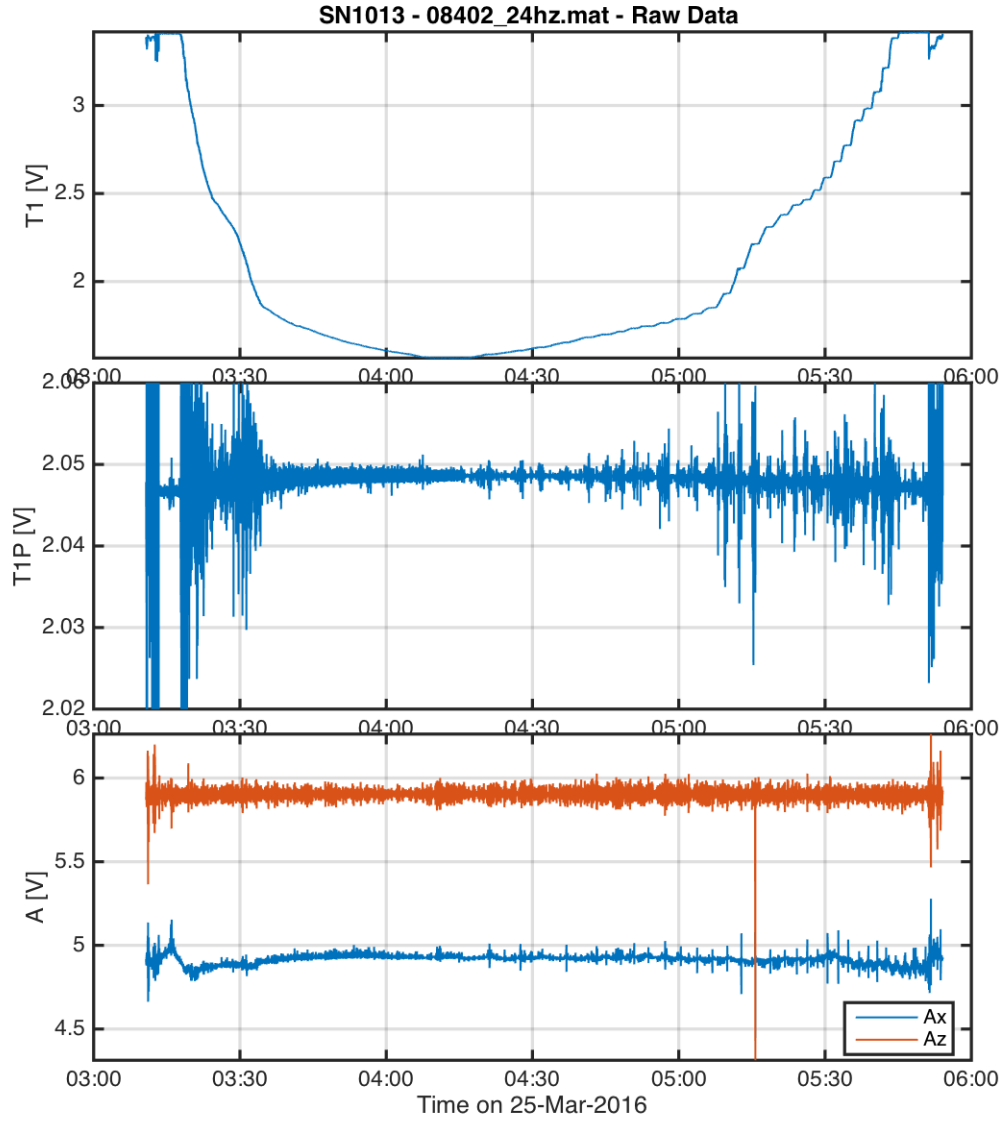


Figure 2: Time-series of raw χ pod data (voltages) for one cast: (a) Temperature (b) dT/dt (c) Accelerometer.

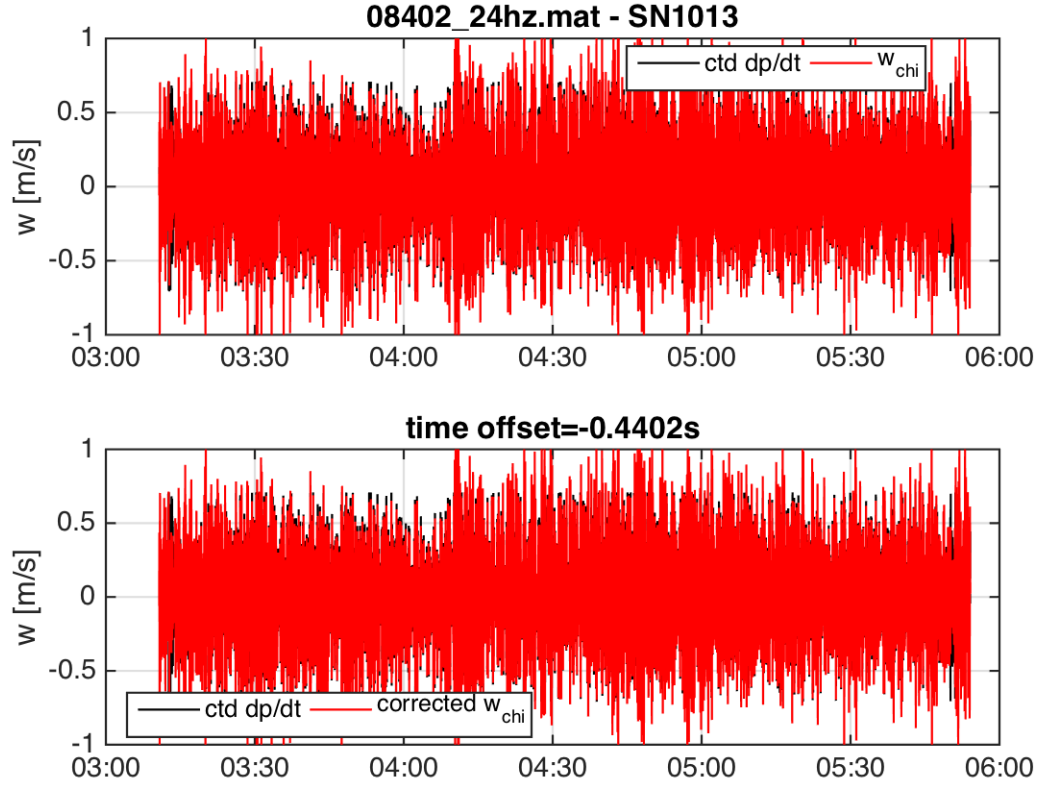


Figure 3: Time series of dp/dt from CTD and integrated vertical accelerations from χ_{pod} . Top is original, bottom is after time-offset applied to χ_{pod} .

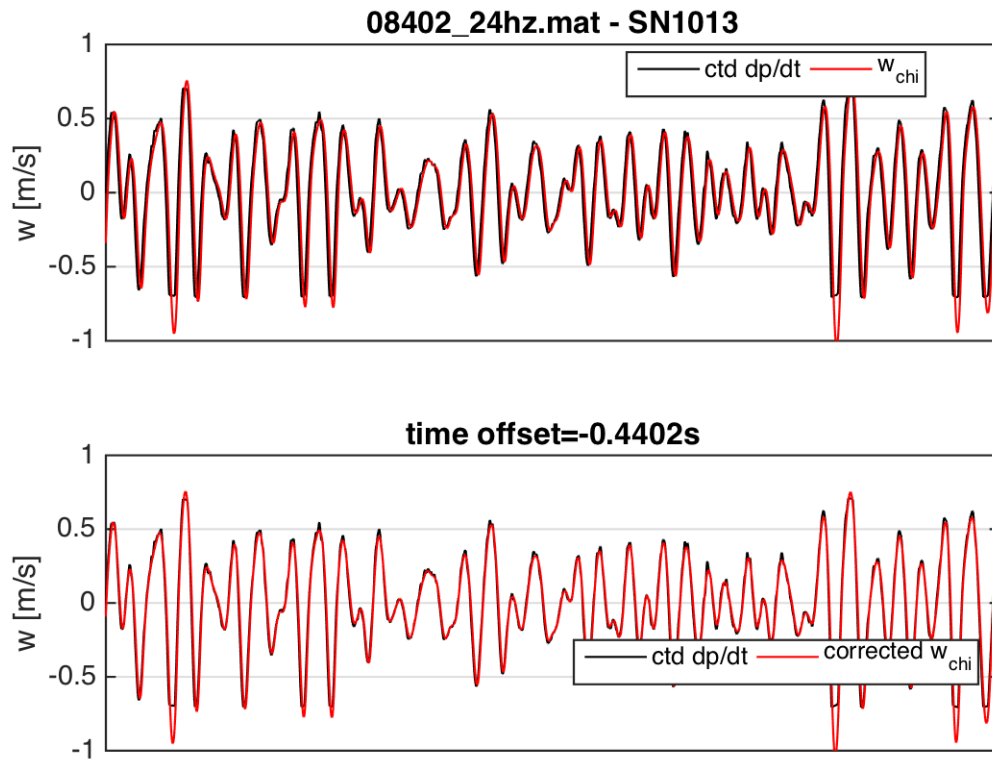


Figure 4: Zoom-in of previous plot allowing one to check if time-offset is correct. Bottom panel should match.

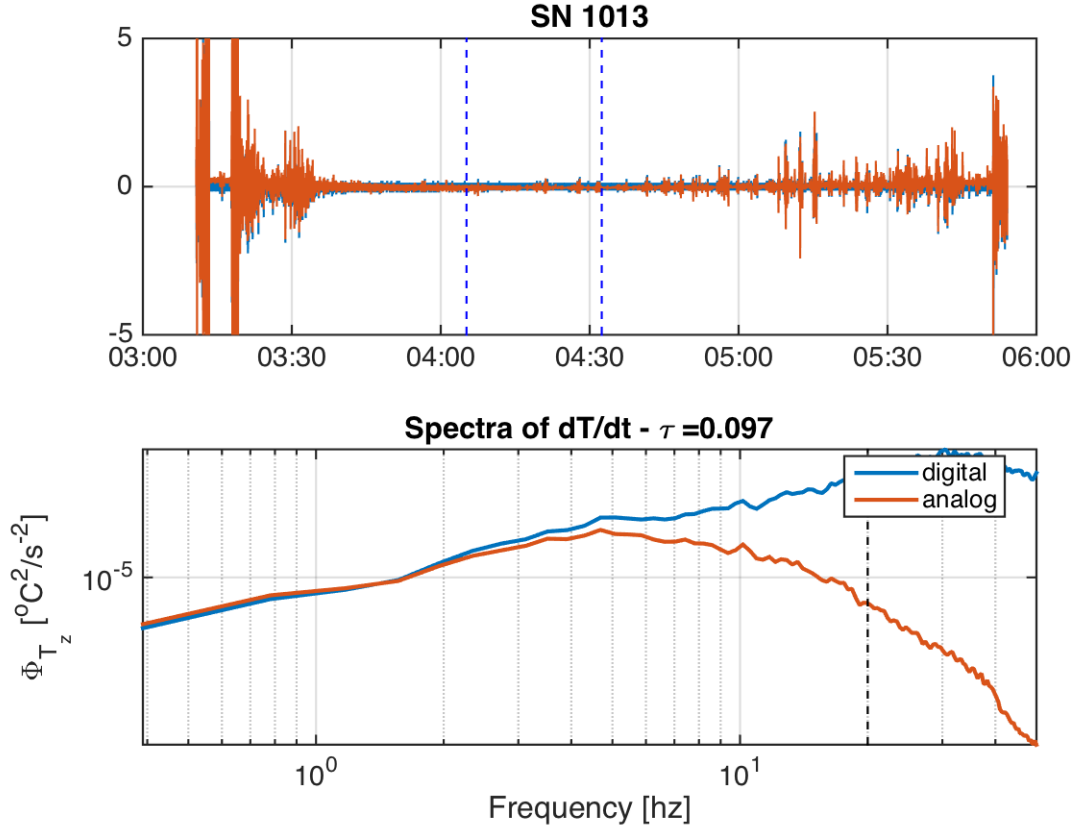


Figure 5: dT/dt spectra for analog (TP) and digital derivatives. Spectra levels should match; if not time-constant may be wrong.

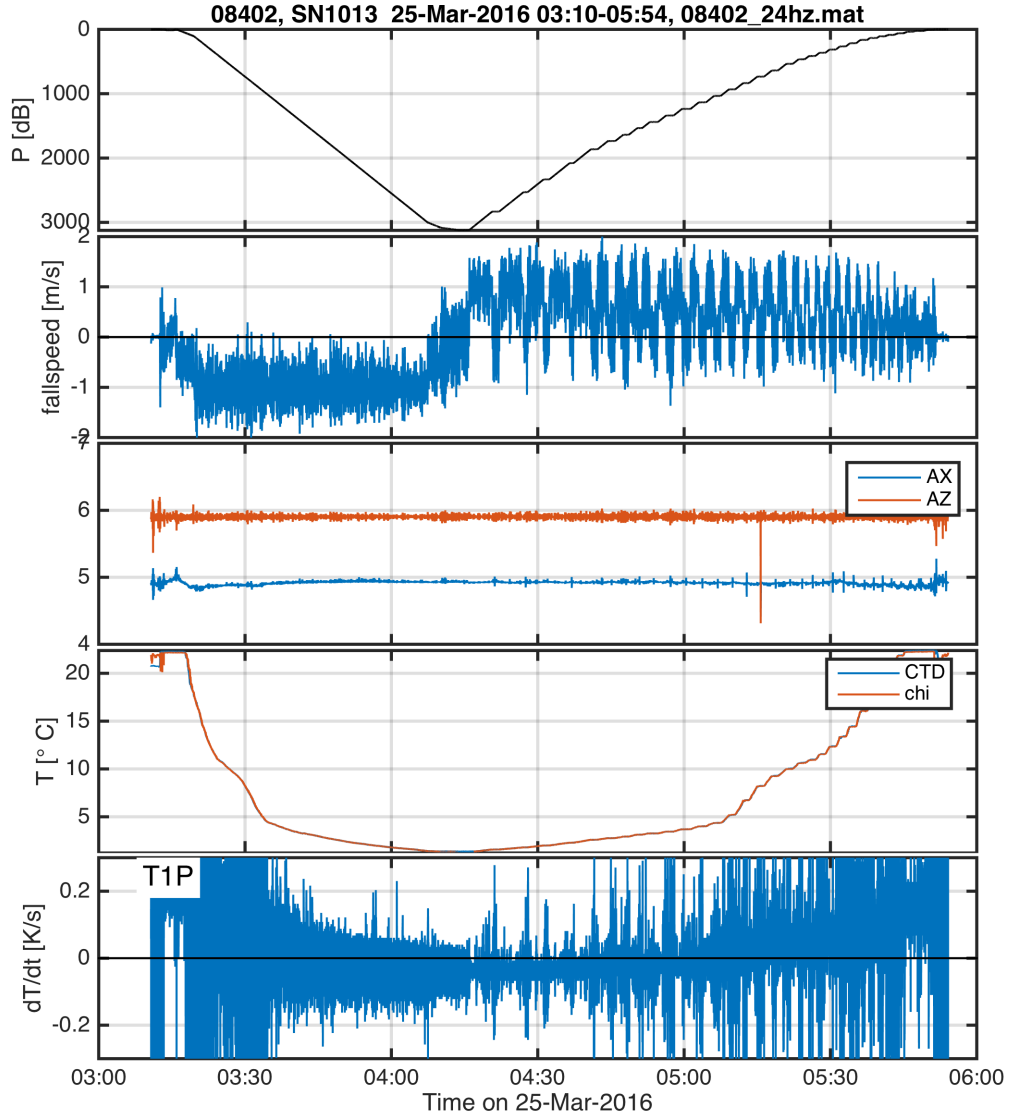


Figure 6: Time series of aligned and calibrated χ pod data for one cast.

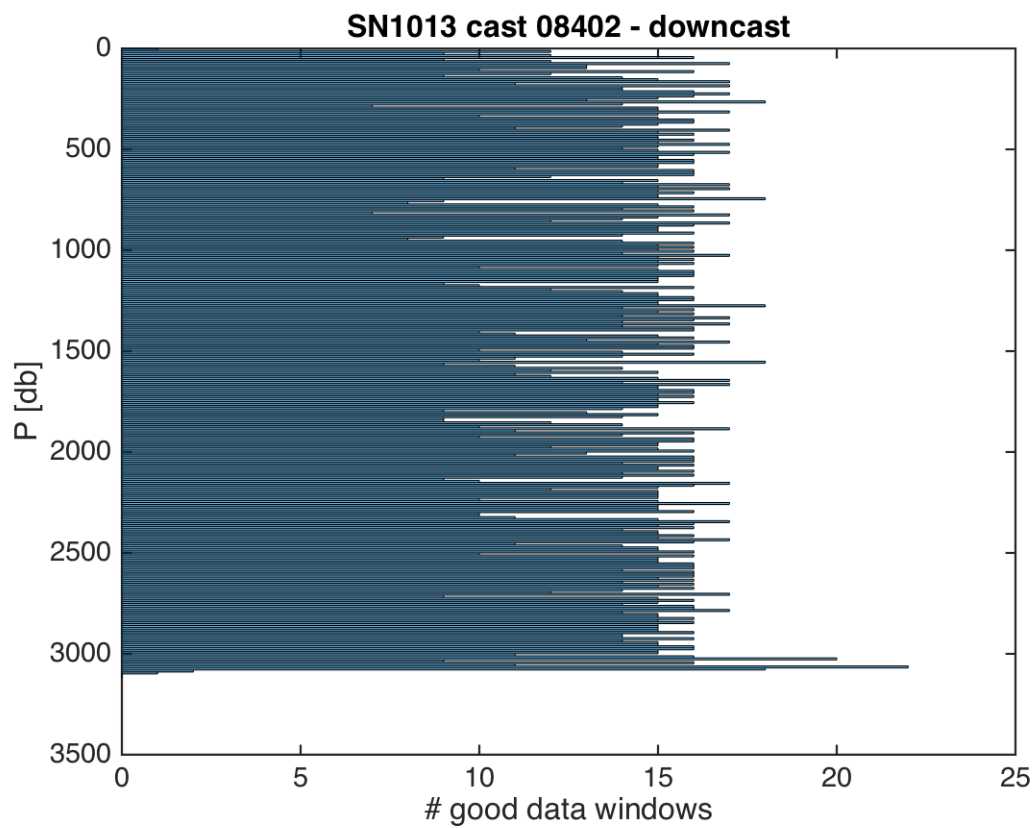


Figure 7: Histogram of the number of good χ pod data windows in each 10m depth bin.

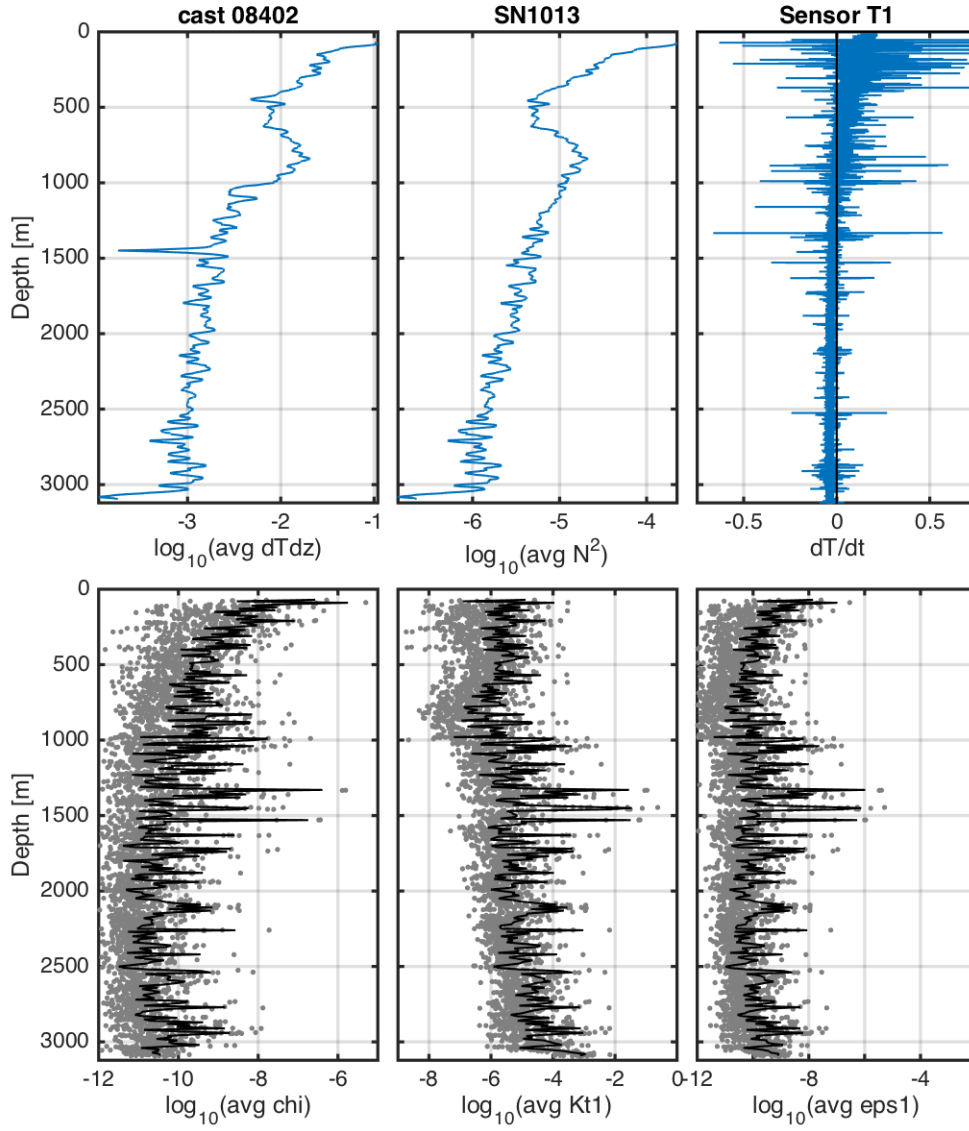


Figure 8: Results of chipod calculation for one cast. Black lines in lower panel are 10m binned averages. Note the spike in dT/dz near 1500m that likely causes false high values of χ .

14 m-files

Processing files are kept in a GitHub repository `/mixingsoftware/CTD_Chipod/mfiles/`

- `AlignChipodCTD` : Align CTD and chipod data.
- `Compute_N2_dTdz_forChi`
- `get_chipod_chi` : Main function that does chipod calculation.
- `load_chipod_data`
- `MakeCtdChiWindows`
-